

SPORTS MEDICINE

PORTFOLIO



ANCHORS
ABLATORS
CANNULAS
ETHOS BUTTON
ETHOS WIRES
ETHOSTAPE
SHAVER BLADES & BURS
INTERFERENCE SCREWS


HNM 
SPORTS MEDICINE

www.hnmmedical.com



TABLE OF CONTENTS

ETHOS LATERAL ANCHOR	3
ETHOS MEDIAL & XT ANCHOR	5
ETHOS PRO LABRAL ANCHOR	7
ETHOS PRO SMALL JOINT ANCHOR	10
ETHOS INTERFERENCE SCREWS	13
DISPOSABLE CANNULA	16
HNM DISPOSABLE RF BIPOLAR	18
ABLATOR	18
ETHOS WIRE	20
ETHOS BUTTON	22
ETHOSTAPE	24
SUTURE PASSER	25
SHAVER BLADES & BURS	26

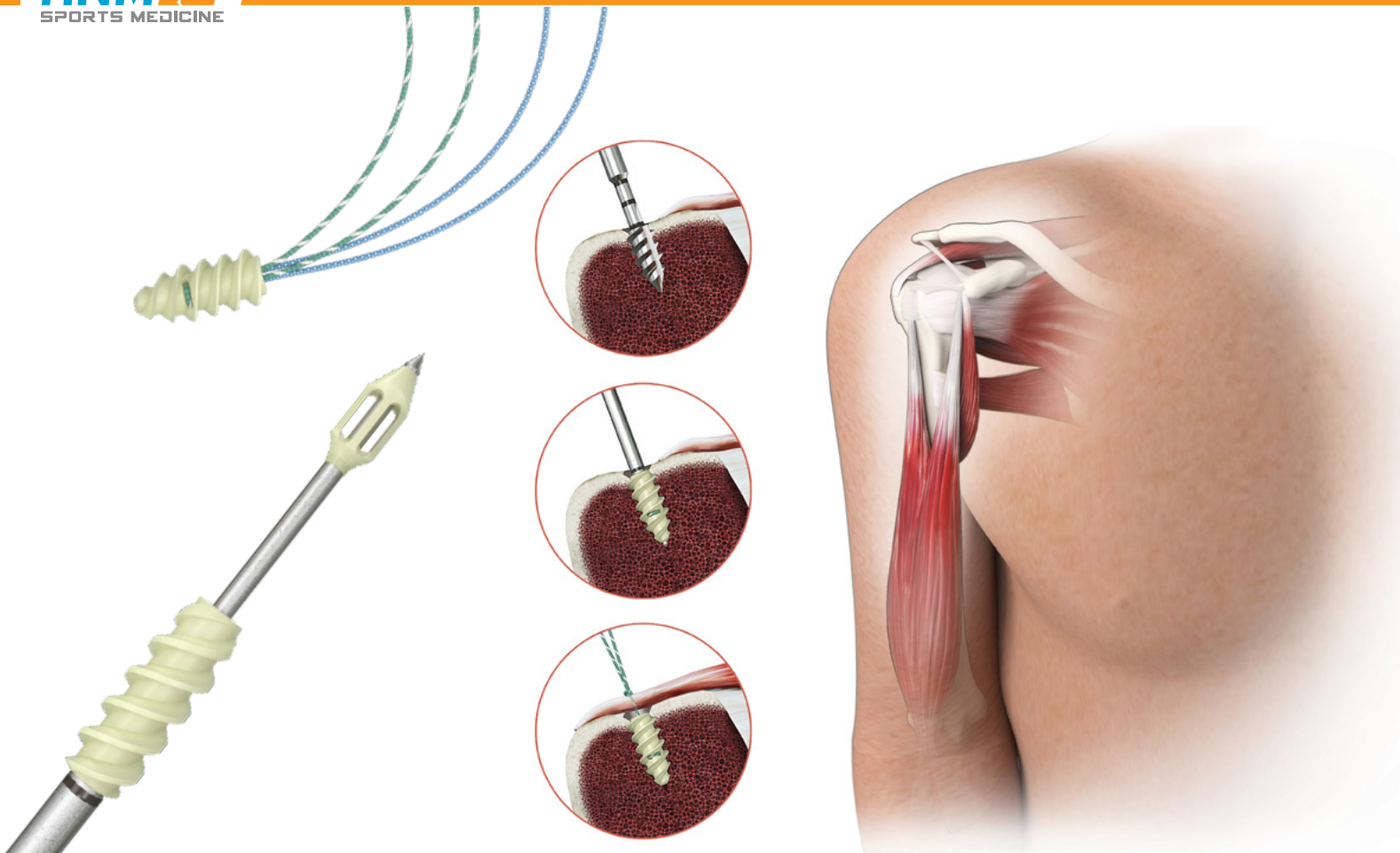

ETHOS LATERAL ANCHOR


PRODUCT INFORMATION

The anchors are used for knotless fixation of ligament, tendon, bone, or soft tissue to bone in knee, shoulder, foot/ankle, elbow, and hand/wrist procedures. These anchors consist of cannulated anchors with integrated suture punch eyelet which leaves no metal in the body if desired. The device is available in varying profile diameters. Each implant is manufactured from Titanium 6AL-4V ELI per ASTM F136 or ZENIVA™ ZA-500 PEEK (polyetheretherketone) from Solvay Advanced Polymers per ASTM F2026.

FEATURES & BENEFITS

- Preloaded fully disposable delivery system with integrated stainless steel punch
- Two piece 'Stable Drive' handle design
- PEEK anchor available in 4.5mm, 5.5mm and 6.5mm diameters
- Dual slotted suture bullet for suture capture from multiple locations and knotless repairs
- Integrated stainless steel punch for insertion with or without predrilling
- Double lead tapered threads for faster delivery and superior pullout strength
- Designed for maximum suture and anchor engagement for superior resistance to pull out.


ORDERING INFORMATION

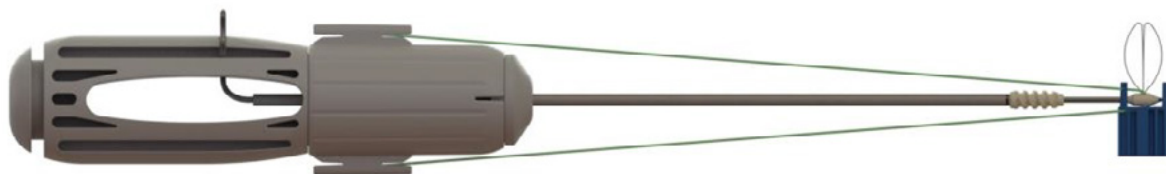

ETHOS LATERAL ANCHORS

ITEM #	DESCRIPTION
HNM-79-4520/PK	ETHOS LATERAL ANCHOR, PEEK, W/PEEK TIP, 4.5mm x 20mm, KNOTLESS, 5/PK
HNM-79-5520/PK	ETHOS LATERAL ANCHOR, PEEK, W/PEEK TIP, 5.5mm x 20mm, KNOTLESS, 5/PK
HNM-79-6520/PK	ETHOS LATERAL ANCHOR, PEEK, W/PEEK TIP, 6.5mm x 20mm, KNOTLESS, 5/PK
HNM-91-4520-TI/PK	ETHOS LATERAL ANCHOR, PEEK, W/TITANIUM TIP, 4.5mm x 20mm, KNOTLESS, 5/PK
HNM-91-5520-TI/PK	ETHOS LATERAL ANCHOR, PEEK, W/TITANIUM TIP, 5.5mm x 20mm, KNOTLESS, 5/PK
HNM-91-6520-TI/PK	ETHOS LATERAL ANCHOR, PEEK, W/TITANIUM TIP, 6.5mm x 20mm, KNOTLESS, 5/PK

INSTRUMENTS

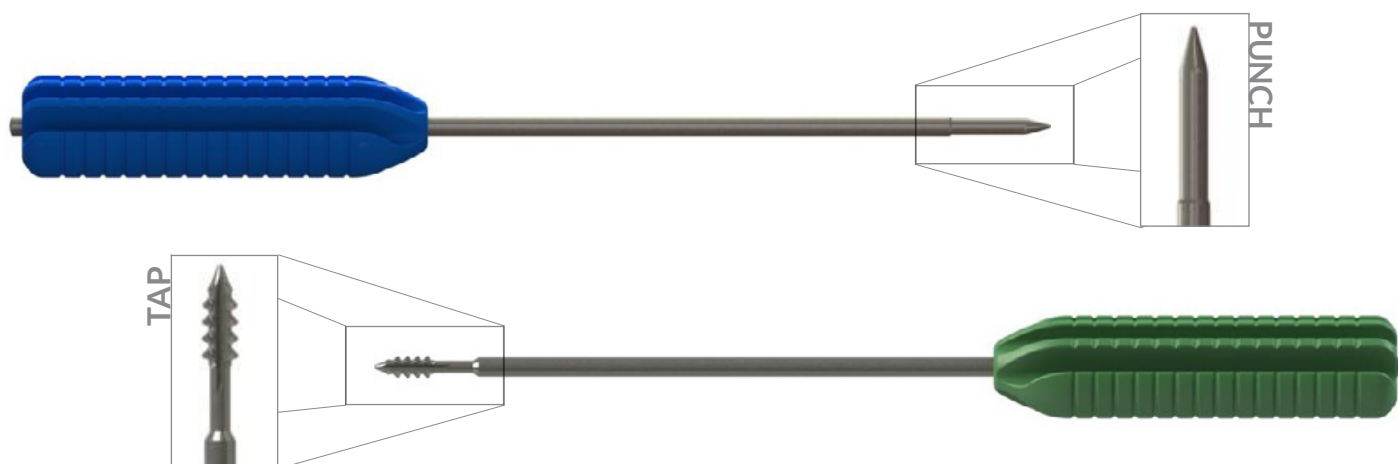

ITEM #	DESCRIPTION
HNM-17-01-45P	ETHOS PUNCH, 4.5mm
HNM-17-01-45T	ETHOS TAP, 4.5mm
HNM-17-01-55P	ETHOS PUNCH, 5.5mm
HNM-17-01-55T	ETHOS TAP, 5.5mm
HNM-17-01-65P	ETHOS PUNCH, 6.5mm
HNM-17-01-65T	ETHOS TAP, 6.5mm



ETHOS MEDIAL & XT ANCHOR


PRODUCT INFORMATION

The anchors are used for fixation of ligament, tendon, bone, or soft tissue to bone in knee, shoulder, foot/ankle, elbow, and hand/wrist procedures. These anchors consist of cannulated anchors with integrated suture attachment. The device is available in varying profile diameters. Each implant is manufactured from Titanium 6AL-4V ELI per ASTM F136 or ZENIVA™ ZA-500 PEEK (polyetheretherketone) from Solvay Advanced Polymers per ASTM F2026.

FEATURES & BENEFITS

- Preloaded fully disposable delivery system with fully contained sutures
- Two piece 'Stable Drive' handle design
- PEEK and Titanium anchor available in 4.5mm, 5.5mm, and 6.5mm diameters
- Integrated stainless steel punch for insertion without predrilling
 - Stainless steel punch removed upon insertion
- Double lead tapered threads for faster delivery
- ETHOS Medial Anchor has two #2 High Strength UHMWPE/Polyester sutures:
 - One Blue suture
 - One White suture or Green suture with White accenture
- ETHOS XT Anchor has three #2 High Strength UHMWPE/Polyester sutures

ORDERING INFORMATION

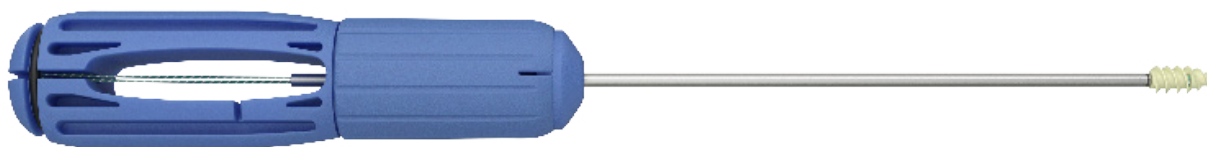


ETHOS MEDIAL ANCHOR

ITEM #	DESCRIPTION
HNM-77-4515/PK	ETHOS MEDIAL ANCHOR, PEEK, 4.5mm x 15mm, (2) #2 SUTURE, 5/PK
HNM-77-5515/PK	ETHOS MEDIAL ANCHOR, PEEK, 5.5mm x 15mm, (2) #2 SUTURE, 5/PK
HNM-77-6515/PK	ETHOS MEDIAL ANCHOR, PEEK, 6.5mm x 15mm, (2) #2 SUTURE, 5/PK
HNM-91-4515-TIPK	ETHOS MEDIAL ANCHOR, TITANIUM, 4.5mm x 15mm, (2) #2 SUTURE, 5/PK
HNM-91-5515-TI/PK	ETHOS MEDIAL ANCHOR, TITANIUM, 5.5mm x 15mm, (2) #2 SUTURE, 5/PK
HNM-91-6515-TI/PK	ETHOS MEDIAL ANCHOR, TITANIUM, 6.5mm x 15mm, (2) #2 SUTURE, 5/PK


ETHOS XT ANCHOR

ITEM #	DESCRIPTION
HNM-78-5515/PK	ETHOS XT ANCHOR, PEEK, 5.5mm x 15mm, (3) #2 SUTURE, 5/PK
HNM-78-6515/PK	ETHOS XT ANCHOR, PEEK, 6.5mm x 15mm, (3) #2 SUTURE, 5/PK
HNM-92-5515-TI/PK	ETHOS XT ANCHOR, TITANIUM, 5.5mm x 15mm, (3) #2 SUTURE, 5/PK
HNM-92-6515-TI/PK	ETHOS XT ANCHOR, TITANIUM, 6.5mm x 15mm, (3) #2 SUTURE, 5/PK

INSTRUMENTS

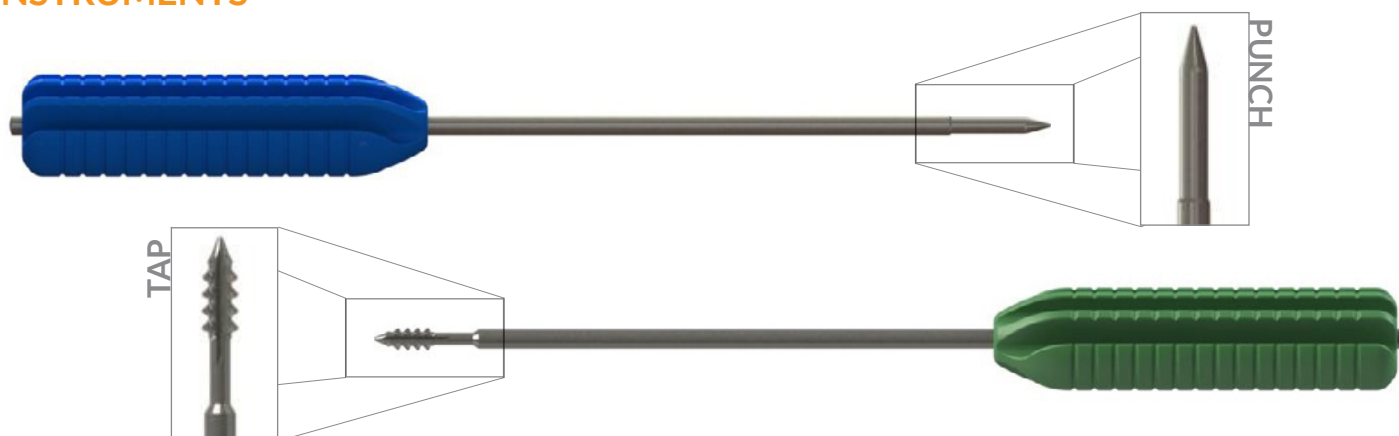

ITEM #	DESCRIPTION
HNM-17-01-45P	ETHOS PUNCH, 4.5mm
HNM-17-01-45T	ETHOS TAP, 4.5mm
HNM-17-01-55P	ETHOS PUNCH, 5.5mm
HNM-17-01-55T	ETHOS TAP, 5.5mm
HNM-17-01-65P	ETHOS PUNCH, 6.5mm
HNM-17-01-65T	ETHOS TAP, 6.5mm

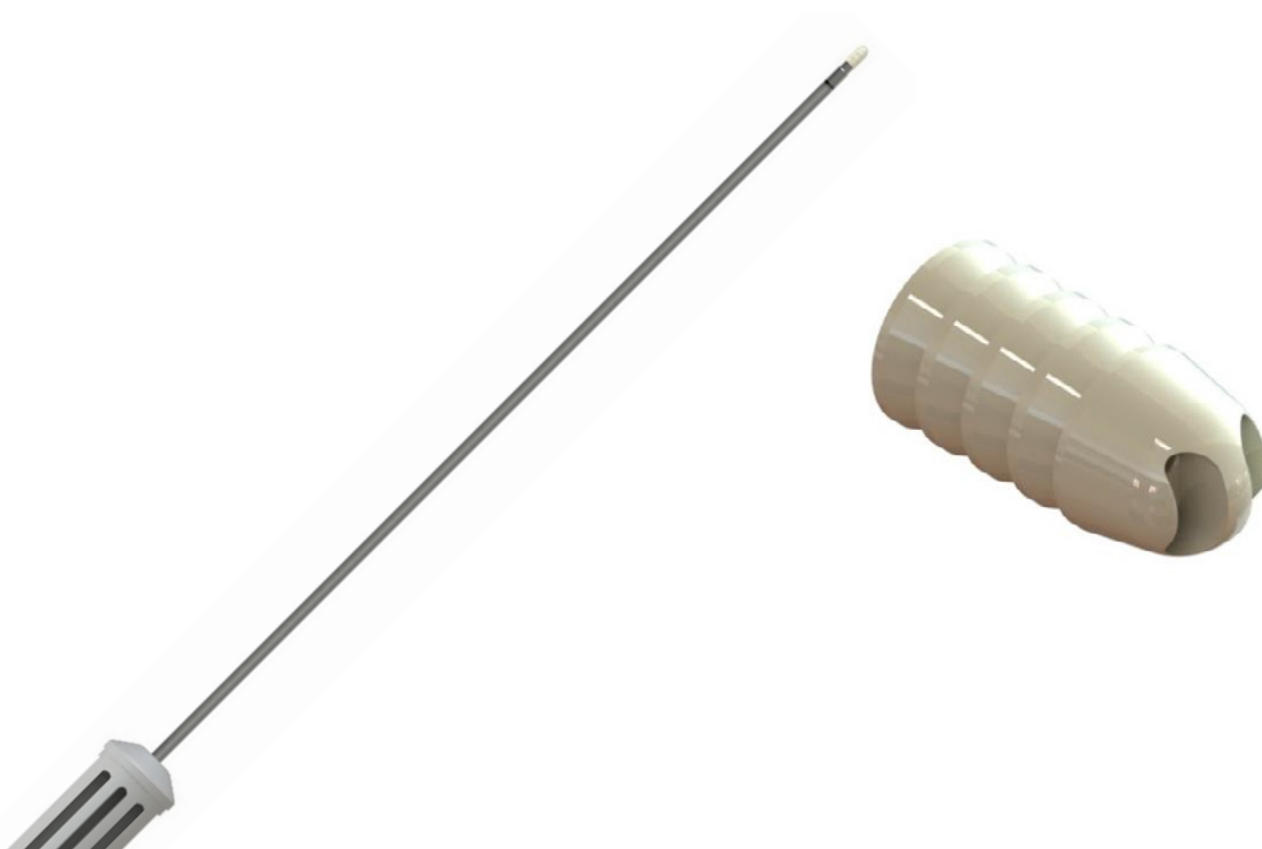

ETHOS PRO LABRAL ANCHOR

PRODUCT INFORMATION

The ETHOS Pro Labral Anchor are for use in fixation of ligament, tendon, bone, or soft tissue to bone in knee, shoulder, foot/ankle, elbow, and hand/wrist procedures. The ETHOS Pro Labral Anchor consists of a single anchor with integrated multiple suture attachments. The anchors are loaded on individual inserters with integrated sutures, sterile, for single use only. The device is available in varying profile diameters and lengths. Each implant is manufactured from Titanium 6AL-4V ELI per ASTM F136 or ZENIVA™ ZA-500 PEEK (polyetheretherketone) from Solvay Advanced Polymers per ASTM F2026.

FEATURES & BENEFITS

- Preloaded fully disposable delivery system with fully contained sutures
- PEEK anchor available in 2.9mm and 3.3mm diameters
- Ribbed body for superior anchor pullout strength
- Instruments designed for ease of delivery
 - No-slip Drill/Delivery Guide for pinpoint delivery and minimal tissue disruption
 - Drills sized for optimal anchor pullout strength
- #2 High Strength UHMWPE/Polyester suture

ORDERING INFORMATION



ETHOS PRO LABRAL ANCHOR


ITEM #	DESCRIPTION
HNM-33-151308/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 2.9mm x 6mm, (2) #2 SUTURE, 5/PK
HNM-33-151164/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 2.9mm x 10mm, (2) #2 SUTURE, 5/PK
HNM-33-151310/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 3.3mm x 6mm, (2) #2 SUTURE, 5/PK
HNM-33-151311/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 3.3mm x 10mm, (2) #2 SUTURE, 5/PK
HNM-33-151312/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 2.9mm x 6mm, (2) #2 SUTURE, 5/PK
HNM-33-151313/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 2.9mm x 10mm, (2) #2 SUTURE, 5/PK
HNM-33-151314/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 3.3mm x 6mm, (2) #2 SUTURE, 5/PK
HNM-33-151315/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 3.3mm x 10mm, (2) #2 SUTURE, 5/PK
HNM-33-151278/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 2.9mm x 6mm, (1) #2 SUTURE, 5/PK
HNM-33-151279/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 2.9mm x 10mm, (1) #2 SUTURE, 5/PK
HNM-33-151280/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 3.3mm x 6mm, (1) #2 SUTURE, 5/PK
HNM-33-151281/PK	ETHOS PRO LABRAL ANCHOR, PEEK, 3.3mm x 10mm, (1) #2 SUTURE, 5/PK
HNM-33-151282/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 2.9mm x 6mm, (1) #2 SUTURE, 5/PK
HNM-33-151283/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 2.9mm x 10mm, (1) #2 SUTURE, 5/PK
HNM-33-151284/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 3.3mm x 6mm, (1) #2 SUTURE, 5/PK
HNM-33-151285/PK	ETHOS PRO LABRAL ANCHOR, TITANIUM, 3.3mm X 10mm, (1) #2 SUTURE, 5/PK

INSTRUMENTS


DRILL WITH DEPTH STOP




ITEM #	DESCRIPTION
HNM-33-151484	ETHOS PRO DRILL W/DEPTH STOP, LABRAL, 2.9mm
HNM-33-151485	ETHOS PRO DRILL W/DEPTH STOP, LABRAL, 3.3mm

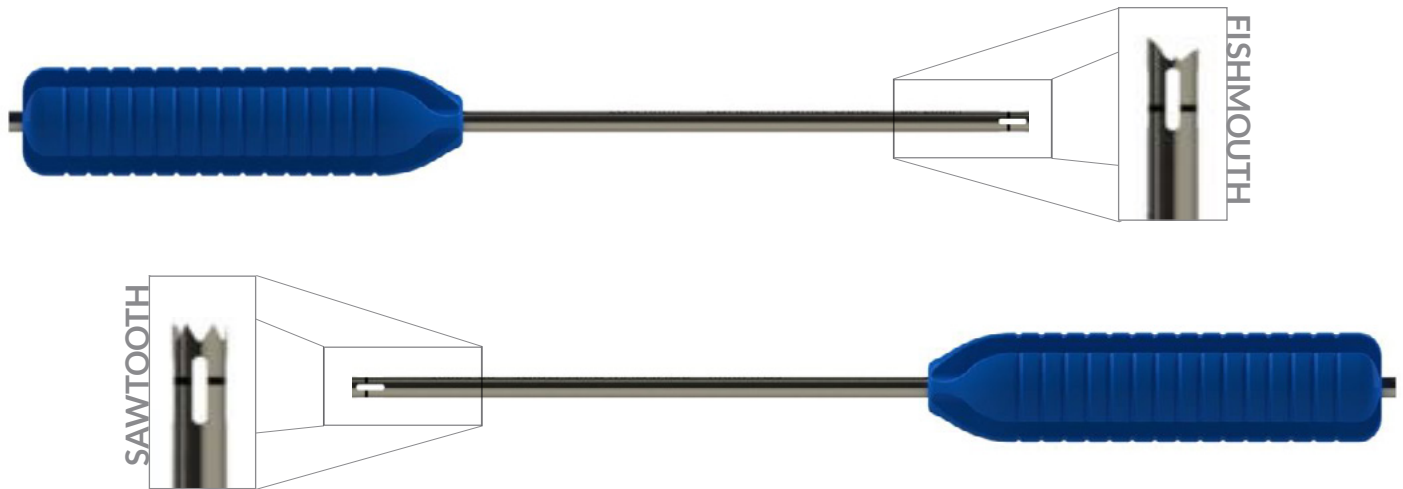

ITEM #	DESCRIPTION
HNM-17-1001-2	ETHOS PRO GUIDE, FISHMOUTH
HNM-33-151489	ETHOS PRO GUIDE, SAWTOOTH

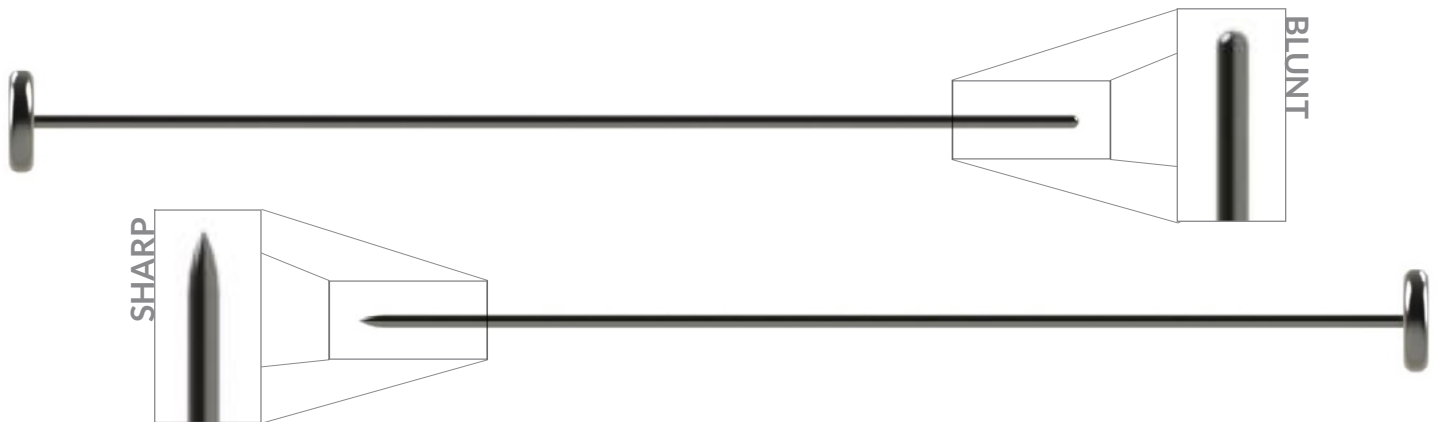

ITEM #	DESCRIPTION
HNM-17-1001-24	ETHOS PRO OBTURATOR, BLUNT
HNM-33-151491	ETHOS PRO OBTURATOR, SHARP

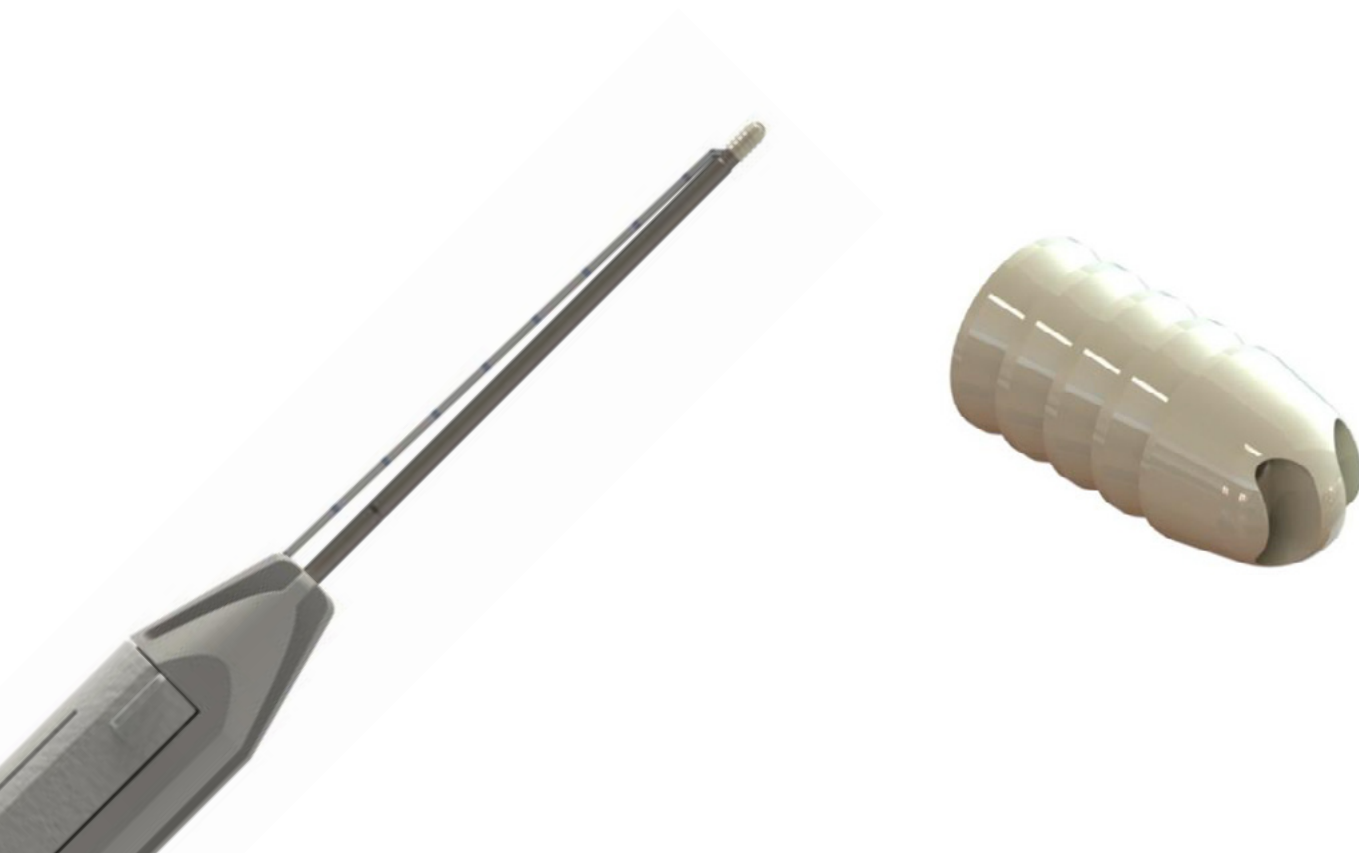

ETHOS PRO SMALL JOINT ANCHOR

PRODUCT INFORMATION

The ETHOS Pro Small Joint Anchor are for use in fixation of ligament, tendon, bone, or soft tissue to bone in knee, shoulder, foot/ankle, elbow, and hand/wrist procedures. The ETHOS Pro Small Joint Anchor consists of a single anchor with integrated multiple suture attachments. The anchors are provided loaded on individual inserters with integrated sutures with needles, sterile, for single use only. The device is available in varying profile diameters and lengths. Each implant is manufactured from Titanium 6AL-4V ELI per ASTM F136 or ZENIVA™ ZA-500 PEEK (polyetheretherketone) from Solvay Advanced Polymers per ASTM F2026.

FEATURES & BENEFITS

- Preloaded fully disposable delivery system with fully contained sutures with needles
- PEEK anchor available in 2.9mm and 3.3mm diameters
- Ribbed body for superior anchor pullout strength
- Shorter shaft for greater stability
- Instruments designed for ease of delivery
 - No-slip Drill/Delivery Guide for pinpoint delivery and minimal tissue disruption
 - Drills sized for optimal anchor pullout strength
- (2) High Strength UHMWPE/Polyester suture and (4) attached needles

ORDERING INFORMATION



ETHOS PRO SMALL JOINT ANCHORS

ITEM #	DESCRIPTION
HNM-33-151443/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 2.9mm x 6mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151444/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 2.9mm x 10mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151445/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 3.3mm x 6mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151446/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 3.3mm x 10mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151451/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 2.9mm x 6mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK
HNM-33-151452/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 2.9mm x 10mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK
HNM-33-151453/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 3.3mm x 6mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK
HNM-33-151454/PK	ETHOS PRO SMALL JOINT ANCHOR, PEEK, 3.3mm x 10mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK
HNM-33-151463/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 2.9mm x 6mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151464/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 2.9mm x 10mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151465/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 3.3mm x 6mm, (2) #2 SUTURE, (4) MO-6 NEEDLE, 5/PK
HNM-33-151466/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 3.3mm x 10mm, (2) #2 SUTURE W/(4) MO-6 NEEDLE, 5/PK
HNM-33-151471/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 2.9mm x 6mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK
HNM-33-151472/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 2.9mm x 10mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK
HNM-33-151473/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 3.3mm x 6mm, (2) #2-0 SUTURE, (4) CT-2 NEEDLE, 5/PK
HNM-33-151474/PK	ETHOS PRO SMALL JOINT ANCHOR, TITANIUM, 3.3mm x 10mm, (2) #2-0 SUTURE W/(4) CT-2 NEEDLE, 5/PK

INSTRUMENTS


DRILL FOR SMALL JOINT



ANGLED GUIDE



ITEM #	DESCRIPTION
HNM-33-151492	ETHOS PRO ANGLED GUIDE
HNM-33-151493	2.9mm ETHOS PRO DRILL, SMALL JOINT
HNM-33-151494	3.3mm ETHOS PRO DRILL, SMALL JOINT

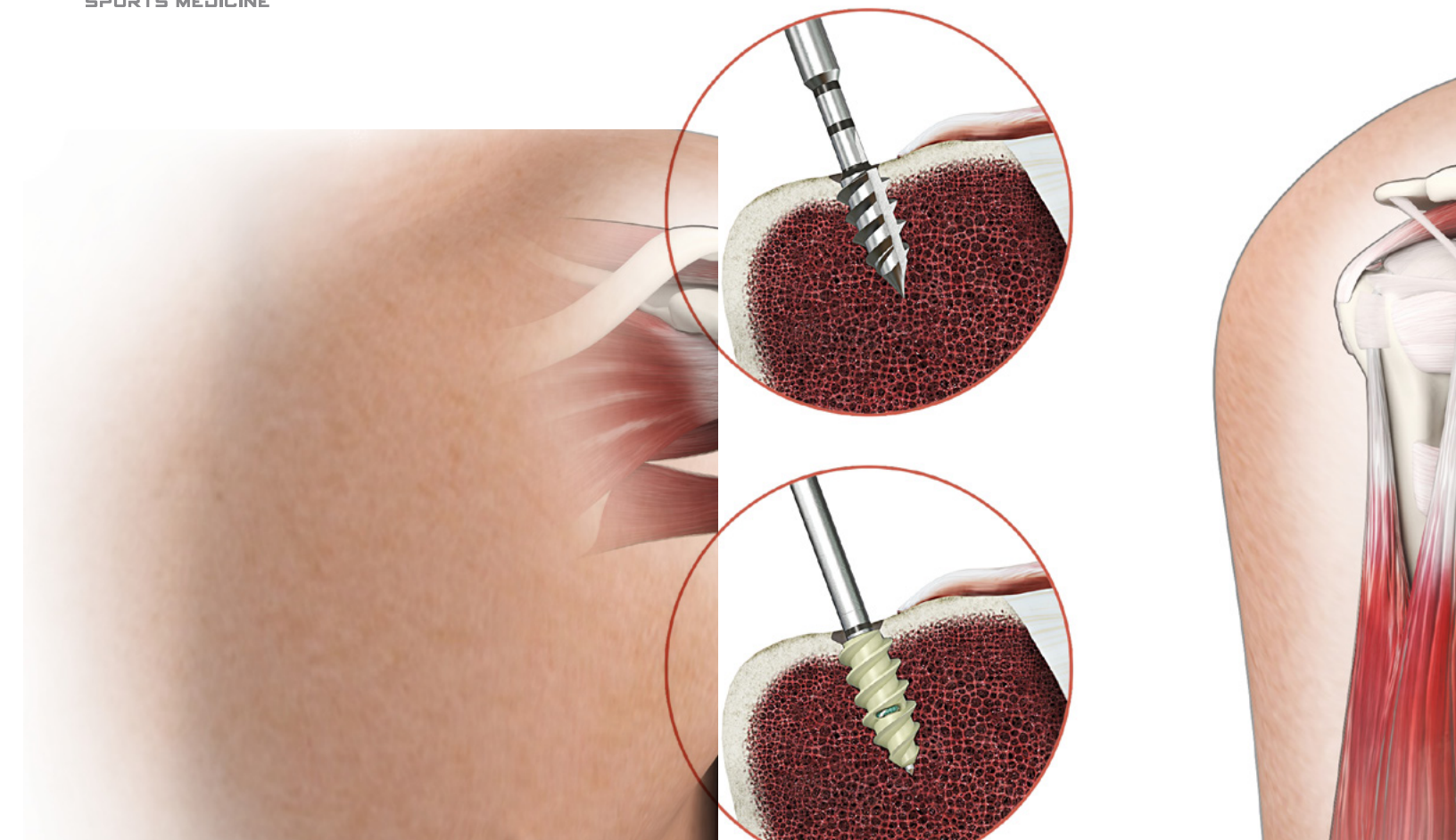

ETHOS INTERFERENCE & MINI INTERFERENCE SCREW

PRODUCT INFORMATION

The ETHOS Interference Screws & Mini Interference Screws are interference screws for use in fixation of ligament, tendon, bone, or soft tissue to bone in knee, shoulder, foot/ankle, elbow, and hand/wrist procedures. The device is available in varying profile diameters and lengths. Each implant is manufactured from Titanium 6AL-4V ELI per ASTM F136 or ZENIVA™ ZA-500 PEEK (polyetheretherketone) from Solvay Advanced Polymers per ASTM F2026.

FEATURES & BENEFITS

- Available in PEEK or Titanium composition fully cannulated
- Double lead tapered threads for faster delivery and superior pullout strength
- Interference Screw available in 7.0mm, 8.0mm, 9.0mm, 10.0mm, 11.0mm, 12.0mm diameters in 23mm and 30mm lengths
 - T20 torx drive feature
- Mini Interference Screw available in 4.5mm, 5.0mm, 5.5mm, 6.0mm, 6.5mm, 7.0mm, 8.0mm diameters
 - 15mm length available for 4.5mm through 7.0mm sizes
 - 12mm length available for 7.0mm and 8.0mm sizes
 - 2mm Hex Drive feature
- Designed for maximum suture/soft tissue and anchor engagement for superior resistance to pull out

ORDERING INFORMATION

ETHOS INTERFERENCE SCREW



ITEM #	DESCRIPTION
HNM-81-0723	ETHOS INTERFERENCE SCREW, PEEK, 7.0mm x 23mm
HNM-81-0730	ETHOS INTERFERENCE SCREW, PEEK, 7.0mm x 30mm
HNM-81-0823	ETHOS INTERFERENCE SCREW, PEEK, 8.0mm x 23mm
HNM-81-0830	ETHOS INTERFERENCE SCREW, PEEK, 8.0mm x 30mm
HNM-81-0923	ETHOS INTERFERENCE SCREW, PEEK, 9.0mm x 23mm
HNM-81-0930	ETHOS INTERFERENCE SCREW, PEEK, 9.0mm x 30mm
HNM-81-1023	ETHOS INTERFERENCE SCREW, PEEK, 10.0mm x 23mm
HNM-81-1030	ETHOS INTERFERENCE SCREW, PEEK, 10.0mm x 30mm
HNM-81-1123	ETHOS INTERFERENCE SCREW, PEEK, 11.0mm x 23mm
HNM-81-1130	ETHOS INTERFERENCE SCREW, PEEK, 11.0mm x 30mm
HNM-81-1223	ETHOS INTERFERENCE SCREW, PEEK, 12.0mm x 23mm
HNM-81-1230	ETHOS INTERFERENCE SCREW, PEEK, 12.0mm x 30mm
HNM-81-1323	ETHOS INTERFERENCE SCREW, PEEK, 13.0mm x 23mm
HNM-81-1330	ETHOS INTERFERENCE SCREW, PEEK, 13.0mm x 30mm

ETHOS MINI INTERFERENCE SCREW



ITEM #	DESCRIPTION
HNM-81-4515	ETHOS MINI INTERFERENCE SCREW, PEEK, 4.5mm x 15mm
HNM-81-5015	ETHOS MINI INTERFERENCE SCREW, PEEK, 5.0mm x 15mm
HNM-81-5515	ETHOS MINI INTERFERENCE SCREW, PEEK, 5.5mm x 15mm
HNM-81-6015	ETHOS MINI INTERFERENCE SCREW, PEEK, 6.0mm x 15mm
HNM-81-6515	ETHOS MINI INTERFERENCE SCREW, PEEK, 6.5mm x 15mm
HNM-81-7012	ETHOS MINI INTERFERENCE SCREW, PEEK, 7.0mm x 12mm
HNM-81-7015	ETHOS MINI INTERFERENCE SCREW, PEEK, 7.0mm x 15mm
HNM-81-8012	ETHOS MINI INTERFERENCE SCREW, PEEK, 8.0mm x 12mm

INSTRUMENTS



ITEM #	DESCRIPTION
HNM-17-0304-040	BONE TUNNEL REAMER, 4.0mm
HNM-17-0304-045	BONE TUNNEL REAMER, 4.5mm
HNM-17-0304-050	BONE TUNNEL REAMER, 5.0mm
HNM-17-0304-055	BONE TUNNEL REAMER, 5.5mm
HNM-17-0304-060	BONE TUNNEL REAMER, 6.0mm
HNM-17-0304-065	BONE TUNNEL REAMER, 6.5mm
HNM-17-0304-070	BONE TUNNEL REAMER, 7.0mm
HNM-17-0304-080	BONE TUNNEL REAMER, 8.0mm
HNM-17-0304-090	BONE TUNNEL REAMER, 9.0mm
HNM-17-0304-100	BONE TUNNEL REAMER, 10.0mm
HNM-17-0304-110	BONE TUNNEL REAMER, 11.0mm
HNM-17-0304-120	BONE TUNNEL REAMER, 12.0mm
HNM-17-0304-130	BONE TUNNEL REAMER, 13.0mm

INTERFERENCE SCREW DRIVER



MINI INTERFERENCE SCREW DRIVER



ITEM #	DESCRIPTION
HNM-18-10003-S	ETHOS INTERFERENCE SCREW DRIVER, 1/4", QUICK CONNECT
HNM-18-10004-L	ETHOS MINI INTERFERENCE SCREW DRIVER, 1/4", QUICK CONNECT, LONG



AXIAL HANDLE

PALM HANDLE



ITEM #	DESCRIPTION
HNM-18-0304-AH	ETHOS QUICK CONNECT AXIAL HANDLE, 1/4"
HNM-18-0304-PH	ETHOS QUICK CONNECT PALM HANDLE, 1/4"



DISPOSABLE CANNULA


PRODUCT INFORMATION

The HNM Disposable Cannula System is a disposable cannula and obturator set for use in general arthroscopic surgical procedures to provide operating portals during the insertion and removal of instruments. The contents are provided sterile, for single use only.

FEATURES & BENEFITS


- Thread depth ensures secure tissue retention in the joint
- Color coded obturators for easy differentiation between sizes and visualization inside the joint (6.5mm - Green, 7.5mm - Blue, 8.5mm - Yellow)
- Large ergonomic obturator handles allow simple insertion
- Transparent Cannula body allows easy instrument and suture visualization
- Double lead fully-threaded body design is gentle on tissue
- Sterile Peel-Packed with 6 per box
- Disposable
- 6.5mm, 7.5mm and 8.5mm inner diameter options
- 70mm and 90mm length options

ORDERING INFORMATION

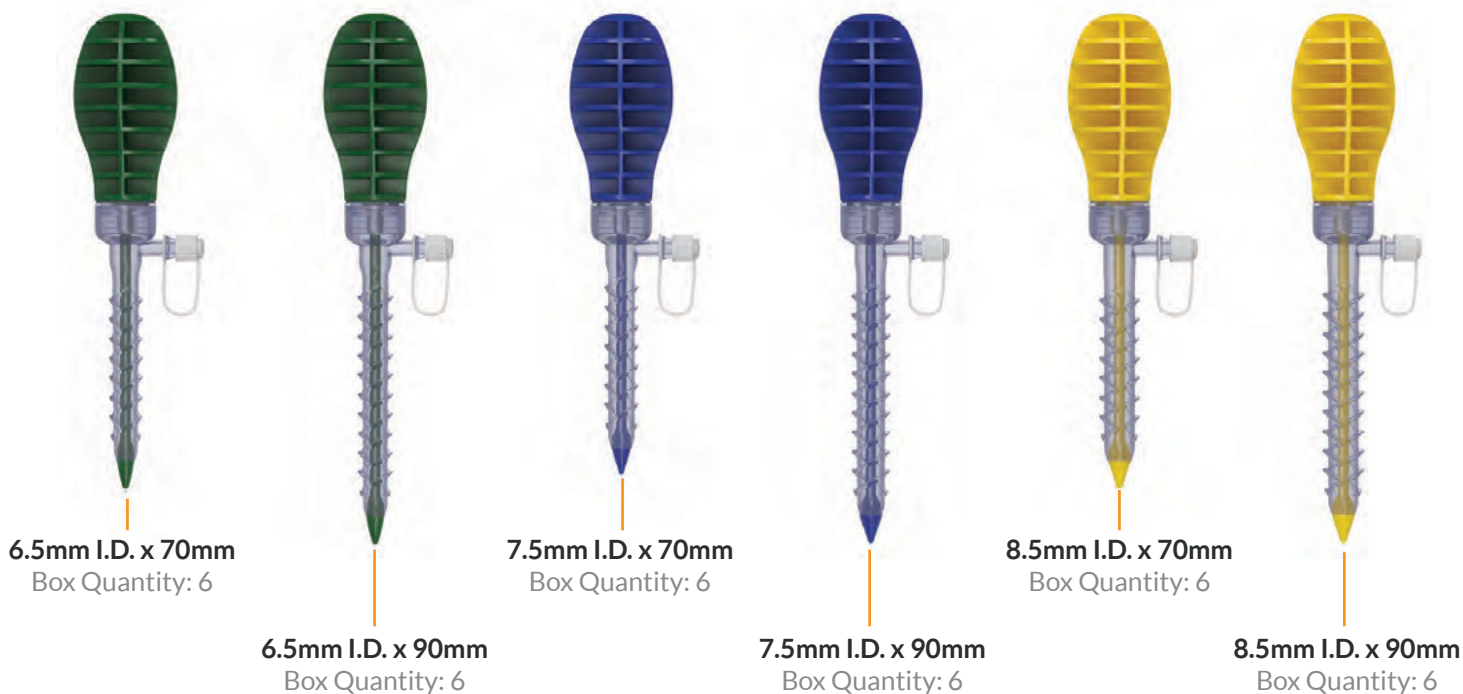

DISPOSABLE CANNULA

ITEM #	DESCRIPTION
HNM-17-100-6570	DISPOSABLE CANNULA, 6.5mm x 70mm, 6/PK
HNM-17-100-6590	DISPOSABLE CANNULA, 6.5mm x 90mm, 6/PK
HNM-17-100-7570	DISPOSABLE CANNULA, 7.5mm x 70mm, 6/PK
HNM-17-100-7590	DISPOSABLE CANNULA, 7.5mm x 90mm, 6/PK
HNM-17-100-8570	DISPOSABLE CANNULA, 8.5mm x 70mm, 6/PK
HNM-17-100-8590	DISPOSABLE CANNULA, 8.5mm x 90mm, 6/PK

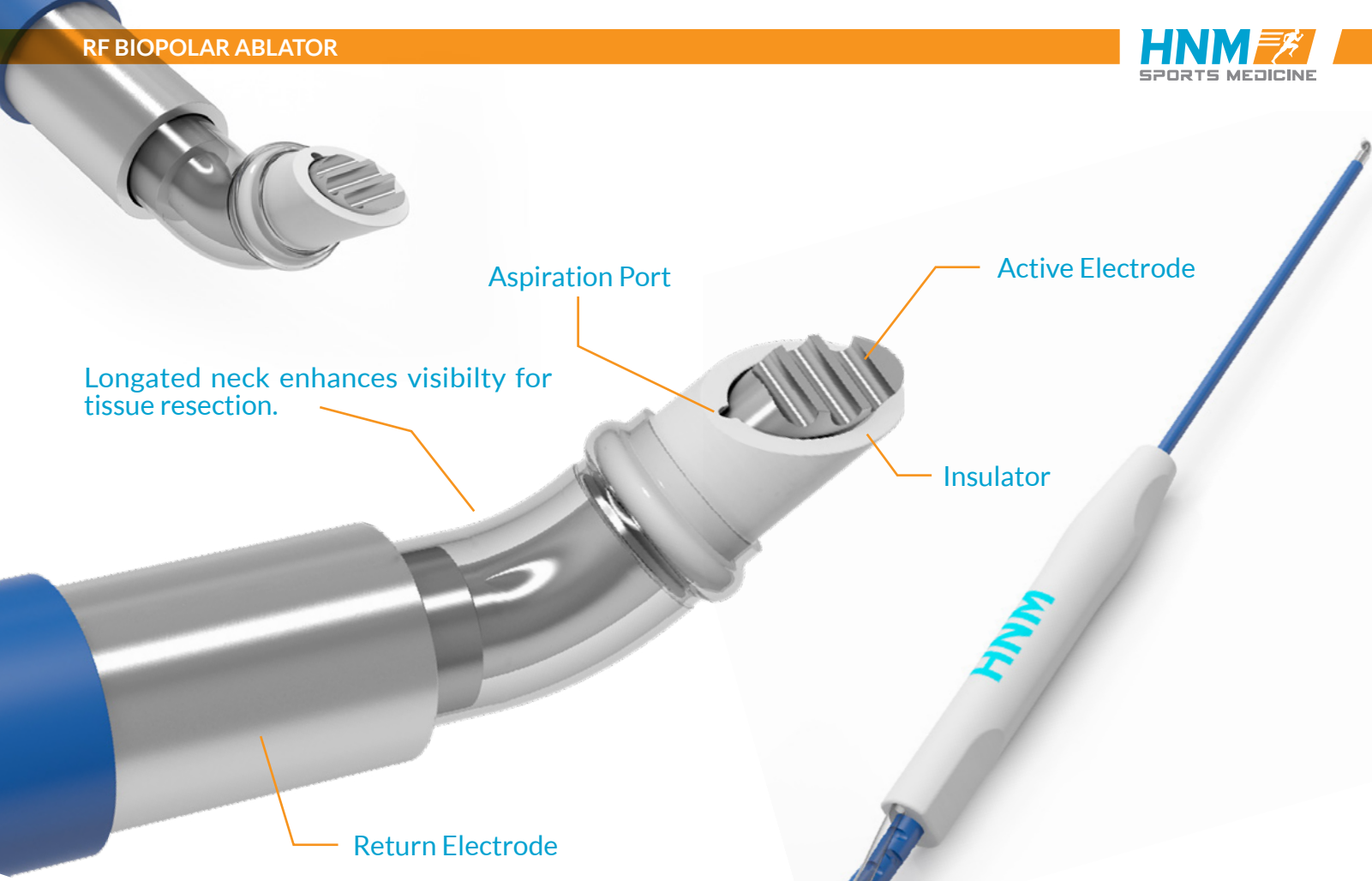

HNM DISPOSABLE RF BIPOLAR ABLATOR

PRODUCT INFORMATION

HNM Disposable Bipolar Ablator is a stainless steel electrode tip, surrounded by a ceramic thermal insulator. The family of bipolar Arthroscopy ablaters by HNM Medical quickly and effectively vaporizes tissue. The family includes sizes for small joint work as well as for shoulders and knees. Aspiration is available on the large size ablaters.

FEATURES & BENEFITS

- Enhanced safety for patients: High efficiency, minimizes risk of thermal injury.
- Sturdier built: Allows for more rugged handling against bone.
- Better aspiration: Aspiration port positioned to prevent clogging.
- Improved visibility: Unique sharp tip design for tissue dissection and angled insulator for tissue elevation.
- Cost efficient: No new equipment needed. Compatible with most standard eletrosurgical generator with footswitch control (i.e. Valleylab, Bovie, Erbe and ConMed).

ORDERING INFORMATION



HNM-31-150365 3.3mm 90° ASPIRATING & NON-ASPIRATING



HNM-31-150367 2.4mm 55° NON-ASPIRATING



HNM-31-150363 3.3mm 50° ASPIRATING & NON -ASPIRATING



HNM-31-150368 1.8mm 60° NON-ASPIRATING

RF BIPOLEAR ABLATOR

ITEM #	DESCRIPTION
HNM-31-150365	3.3mm 90° ASPIRATING & NON-ASPIRATING
HNM-31-150363	3.3mm 50° ASPIRATING & NON -ASPIRATING



ETHOS WIRE

PRODUCT INFORMATION

ETHOS WIRE is cleared by the FDA for use in general soft tissue approximation and/or ligation, including use in cardiovascular procedures, and the use of allograft tissue for orthopedic procedures.

Size 5 through 4-0 (USP) available on spools, pre-cut lengths, needle or non-needle, sterile or non-sterile packs.

FEATURES & BENEFITS

- Four enhanced colors with triaxial pattern design. Differentiates braid for maximum visibility
- 10x stronger than steel on a weighted basis.
- New eBraid design is supple, soft for better handling.
- New core-loc design means stronger knots.

ORDERING INFORMATION

NEEDLE	COLOR	LENGTH	SIZE 25
MO-6 26mm	White Blue	36"-90cm	HNM-e-2C-NN
1/2 Circle			HNM-e-2C



ETHOS BUTTON



PRODUCT INFORMATION

The ETHOS Button CL; a product that is projected to have comparable strength to similar devices indicated for anterior cruciate ligament (ACL) femoral cortical fixation repair.

FEATURES & BENEFITS

- The ETHOS Button CL consists of an ultra-high molecular weight polyethylene (UHMWPE) loop, combined with a titanium plate (button)
- ETHOS Button CL is an suture polyethylene suture offered in a loop configuration. It is truly a champion in strength and durability
- Available in varies sizes from 15mm to 60mm and a adjustable loop

ORDERING INFORMATION

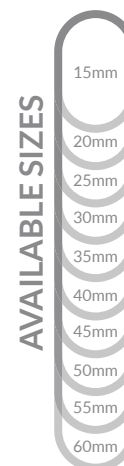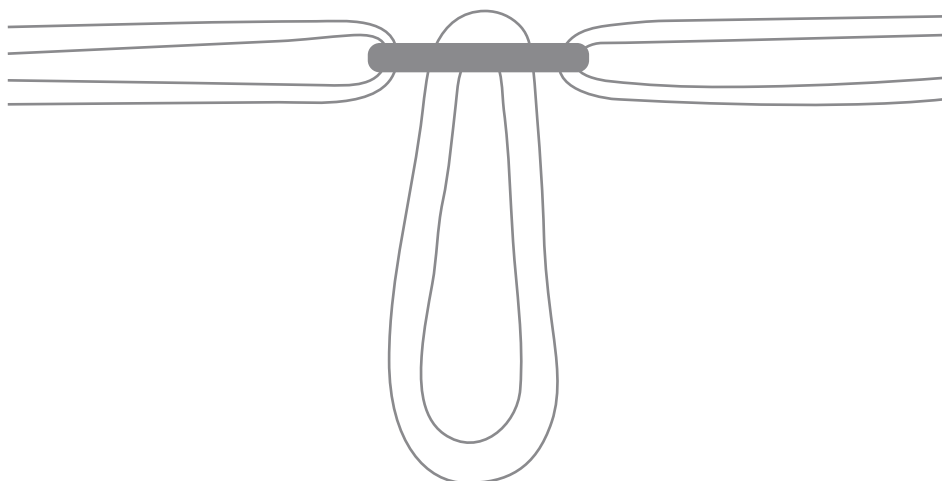

ETHOS BUTTON CL

ITEM #	DESCRIPTION
HNM-OBCL-TI-15	ETHOS BUTTON CL, 15mm LOOP, (2) #5 SUTURE
HNM-OBCL-TI-20	ETHOS BUTTON CL, 20mm LOOP, (2) #5 SUTURE
HNM-OBCL-TI-25	ETHOS BUTTON CL, 25mm LOOP, (2) #5 SUTURE
HNM-OBCL-TI-30	ETHOS BUTTON CL, 30mm LOOP, (2) #5 SUTURE



ETHOS BUTTON AL

ITEM #	DESCRIPTION
HNM-OBAL-80-TI	ETHOS BUTTON AL (ADJUSTABLE LOOP), (2) #5 SUTURE



ETHOSTAPE

PRODUCT INFORMATION


ETHOSTAPE is a Ultra High Molecular Weight (UHMWPE) polyethylene non-absorbable suture offered in a tape configuration.

FEATURES & BENEFITS


- The flat surface of the tape helps brace and distribute more consistent loads while still giving the strength needed
- The broad, flat, braided portion also provides for greater surface area fixation compared to traditional suture
- The colored warp strand enhances visibility
- Our ETHOSTAPE products are great for: Rotator Cuff Repair and Anchor Type Devices

ORDERING INFORMATION

ITEM #	COLOR	SUTURE SIZE	NEEDLE	LENGTH
HNH-209	White / Blue	2mm	Non-Needled	39" (99cm)
HNH-211	White / Black	2mm	Non-Needled	39" (99cm)

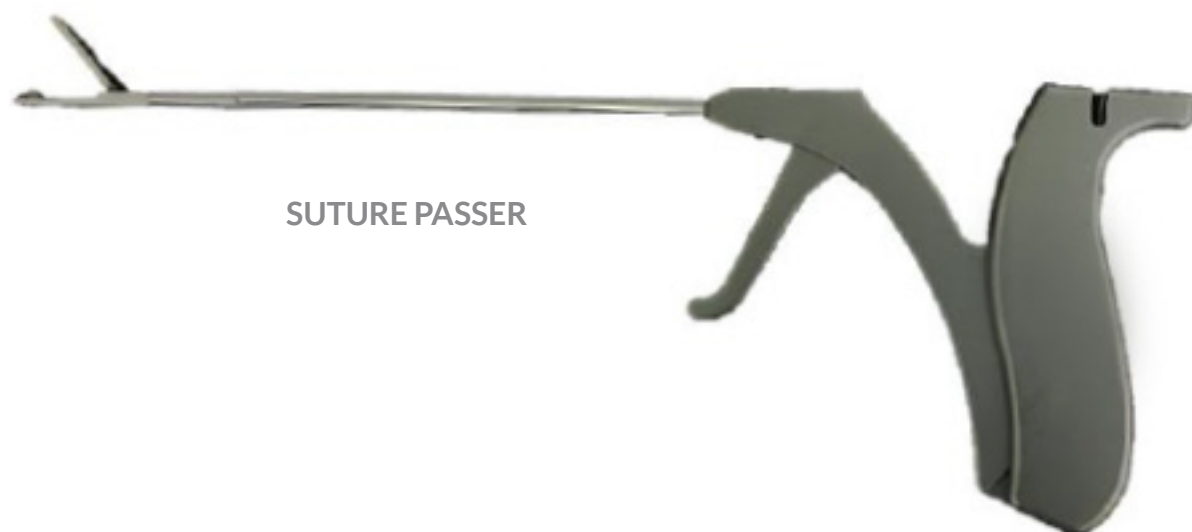

SUTURE PASSER

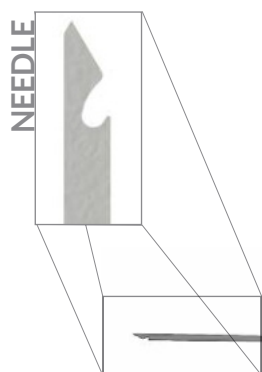

NEEDLE

SUTURE NEEDLE

SUTURE PASSER

PRODUCT INFORMATION

Suture Passer is made from Stainless Steel with side-loading. Suture Needle is made of Nitinol. When needle is deployed, it extends to 1.5cm while other suture passers extend to 1cm.

ORDERING INFORMATION

ITEM #	DESCRIPTION
HNM-33-151344	Suture Passer
HNM-49-151385	Disposable Needle for Suture Passer, Nitinol, Sterile, 10/pk



SHAVER BLADES & BURS

PRODUCT INFORMATION

Disposable Arthroscopic Shaver Blades & Burs for use with Smith & Nephew Dyonics®, ConMed Linvatec® and Stryker® Shaver Systems. A NEW blade for use with your current shaver system.

All blades and burs feature high performance, proprietary plating for reduced friction, strength and precision sharpness.

FEATURES & BENEFITS

- Providing optimal performance with a value strategy.
- No new capital equipment. Just connect and resect.
- New precision blades, sharp every time, clean every time. No proteins, no damage, no doubts.
- Meet cost containment pressures with a product strategy designed to meet your financial objectives.
- Sizes and styles for the most common knee and shoulder procedures.
- No change to inventory practices. Six blades per box.

SHAVER BLADES & BURS COMPATIBLE WITH:

Smith & Nephew Dyonics Handpieces: EP-1, Dyonics Powermax, Dyonics Powermac Elite

Ref. Catalog No's: 7205355, 7205354, 7210542, 7210586



Conmed Linvatec Handpieces: Ergo, Advantage, Advantage Turbo, Apex


Ref. Catalog No's: D4240, D4200, D9824, D9920, C9820, C9824, C9828

Stryker Handpieces: Formula

Ref. Catalog No's: 375-708-500, 375-704-500, 375-701-500

ORDERING INFORMATION

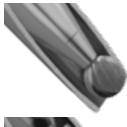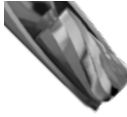BLADE & BUR		COMPARE TO SMITH & NEPHEW DYONICS
	FULL RADIUS BLADE	Full Radius, Full Radius Bonecutter, Synovator or TurboTrimmer Blades
	3.5mm, HNM-32-150370	7205305, 7210486
	4.5 mm, HNM-32-150372	7205306, 7205315, 7206011, 7206013
	5.5 mm, HNM-32-150374	7205307, 7206010, 7206012
	RAPIDCUT BLADE	Incisor, FR Elite, Incisor Plus Elite or Chondroplasty Blades
	3.5mm, HNM-32-150371	7209211, 7205312, 72200095
	4.5 mm, HNM-32-150373	7205313, 7205345, 7209489, 7210499, 7210976
	5.5 mm, HNM-32-150375	7205314, 72200081
	NEATCUT BLADE	TurboTrimmer, Incisor, Incisor Plus, Full Radius Plus, Full Radius Elite, or Trimmer Blades
	3.5mm, HNM-32-150385	7210751
	4.5 mm, HNM-32-150380	7205315, 7209201, 7209202, 7209489, 7210499
	5.5 mm, HNM-32-150386	7205459, 7210909
	ROUND BUR	Notchblaster, Abrader or Shielded Burs
	4.5 mm, HNM-32-150376	7205324, 7205328, 72200722, 72200726
	5.5 mm, HNM-32-150378	7205325, 7205329, 72200723, 72200727
	OVAL BUR	Stonecutter, Acromionizer, Acromioblaster or Tapered Burs
	4.5 mm, HNM-32-150377	6900988, 7205326, 7205330, 7205668, 7209512, 72200724, 72200728, 72200730
	5.5 mm, HNM-32-150379	6900996, 7205327, 7205331, 7205669, 72200725, 72200729, 72200731

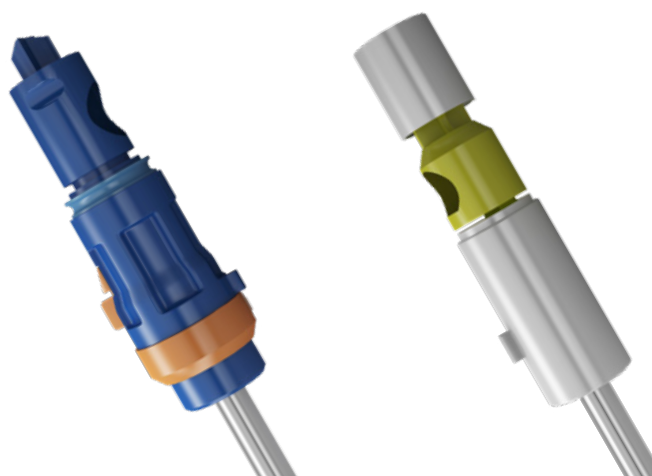

BLADE & BUR		COMPARE TO LINVATEC
	FULL RADIUS BLADE	Full Radius Resector, UltraCut, Tiger or Meniscus Cutter Blades
	3.5mm, HNM-32-150387	35ULT-ST-ZZ, C9248, C9240, C9760
	4.5 mm, HNM-32-150389	9247A, C9405A, C9245, C9242, 9761A, C9405
	5.5 mm, HNM-32-150391	9246A, C9415, C9762
	RAPIDCUT BLADE	Cuda, Mako or Great White Blades
	3.5mm, HNM-32-150388	9399A, C9253, 9297A, 9298A, 9299A, C9254, C9258
	4.5 mm, HNM-32-150390	9599A, C9255, C9298A
	5.5 mm, HNM-32-150392	
	NEATCUT BLADE	Gator or Meniscus Cutter Blades
	3.5mm, HNM-32-150402	C9264, C9760, C9961
	4.5 mm, HNM-32-150397	C9262, 9263A, 9761A
	5.5 mm, HNM-32-150403	9260A, C9762
	ROUND BUR	Cyclone, Oval or Spherical Burs
	4.5 mm, HNM-32-150393	H9111, H9118
	5.5 mm, HNM-32-150395	H9112, H9112RH, H9119
	OVAL BUR	Vortex, Tapered, Oval or Lightning Burs
	4.5 mm, HNM-32-150394	H9101, H9121, H9131, H9101RH, H9121RH, H9133
	5.5 mm, HNM-32-150396	H9102, H9122, H9132, H9103RH

BLADE & BUR		COMPARE TO STRYKER
	FULL RADIUS BLADE	Formula Resector or Scalloped Blades
	3.5mm, HNM-32-150404	375-532-000
	4.5 mm, HNM-32-150406	375-542-000, 375-546-000, 375-552-000
	5.5 mm, HNM-32-150408	375-562-000, 375-566-000
	RAPIDCUT BLADE	Formula Tomcat Blades
	3.5mm, HNM-32-150405	375-535-000
	4.5 mm, HNM-32-150407	375-545-000, 375-555-000
	5.5 mm, HNM-32-150409	375-565-000
	NEATCUT BLADE	Aggressive Plus or Scalloped Blades
	3.5mm, HNM-32-150414	375-534-000
	4.5 mm, HNM-32-150415	375-554-000, 375-544-000, 575-546-000
	5.5 mm, HNM-32-150420	375-564-000, 375-566-000
	ROUND BUR	Formula Round Burs
	4.5 mm, HNM-32-150410	375-940-000, 375-940-012
	5.5 mm, HNM-32-150412	375-950-000, 375-950-012, 375-950-100, 375-950-112
	OVAL BUR	Formula Barrel, Auger or Egg Burs
	4.5 mm, HNM-32-150411	375-941-000, 375-941-012
	5.5 mm, HNM-32-150413	375-450-500, 375-951-000, 375-951-012, 375-951-100, 375-951-112, 375-952-000

NOTES

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

NOTES

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

NOTES

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

NOTES

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

